




MEDIPLUS
INDIA LIMITED

LUER LOCK

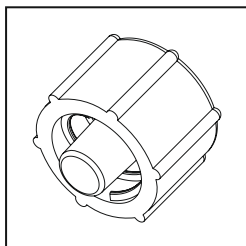

SRN. No: **IN-MF-000007279**
Basic UDI-DI #: **890602530LCZR**

INSTRUCTIONS FOR USE

Materials Used: HDPE

Device General Description: It is a single end luer lock which is used as a male stopper.

Intended Purpose/use: The Sterile luer lock provides a non-injectable end cap. Used to cap off female luer. It is not having access port. The duration of use of this product is 72 hours.

Intended Users: Uses is restricted to a qualified Doctor or a Paramedic.

Principle of operation/Mode of Action: • Sterile Luer Lock is a single end luer lock which is used as a male stopper.

Residual Risk: There is no residual risk associated with the device, however preventive measures are considered in contraindication.

Instructions For Use: • Inspect visually to ascertain that package is intact. • Remove Luer Lock from the sterile packing.

• Connect the device to serve as protective stopper for any standard female luer lock port. • Fix the Luer Lock by rotating it.

Indication: • It is a male luer lock. • Ideal for use with any standard female luer port. • The device provides a non-injectable end.

Contraindications: • Product should not be used in patients with known hypersensitivity to any of the materials used. • Not to be used with NON-STANDARD female luer/mating parts.

Warnings: • The use of this product is restricted to a qualified doctor or a paramedic. • Read Instruction before use. • The product should be used according to the instructions for use. • The product should not be reprocessed. • Visually inspect and carefully check the damage of the device or packaging before use. • Improper transportation and handling may cause structural and/or functional damage to device or packaging. • The product is guaranteed non-toxic, sterile &

non-pyrogenic if the package has not been opened or damaged. • Do not clean or resterilise. • For single use only. • Discard after single use. • If the device is reused, than it may cause bacterial infection to the patient. • Store in cool & dry place. • Do not expose to heat or direct sunlight. • The product should be used immediately after opening the packaging.

Clinical Risk: None

Caution: DISPOSE OFF SAFELY AFTER USE AWAY FROM LIVING BEING.

Note: If the above said caution are not followed by the Paramedic Staff, it can lead to halt the intended clinical benefits to patient.

Patient Population: All age group.

Safer Mode of Use: Up to 72 hrs. for best performance.

Clinical Benefits: Alleviation of disease by providing fluids/medicaments to the patients with channeling other compatible medical device.

Storage: Store in cool & dry place. Do not expose to heat or direct sunlight.

Incidents Report: In case of incidence Immediately inform to Manufacturer (info@mediplusindia.com) and Competent authority of country.

MEDIPLUS (INDIA) LIMITED DISCLAIM ANY RESPONSIBILITY FOR POSSIBLE CONSEQUENCES RESULTING FROM IMPROPER USE.

REF	Catalogue number	Mfg. date with country of manufacture	Manufacturer
Use-by date		Do not use if package is damaged	Caution
LOT	Batch code	STERILE EO Sterilized using ethylene oxide	Do not re-use
Do not re-sterilize		Consult instructions for use	Non-pyrogenic
Single sterile barrier system		Product is not made with natural rubber latex	Does not contain DEHP
MD	Medical Device	UDI Unique Device Identifier	Distributor
Importer		Keep away from sunlight	Keep dry
EU REP	Authorised Representative in the European Union	CE mark European Conformity	NI/IFUS_01, Rev. 07 (Issued date: 16/05/2026)



MEDIPLUS
INDIA LIMITED

CE 1639

Mfd. by:- Mediplus (India) Limited

Unit - 1:- 1261-1262, M.I.E., Part-B, Bahadurgarh - 124507, Haryana (INDIA). M.L.No.: MFG/MD/2018/000068
Unit - 2:- 32, Sector - 16, HSIDC, Bahadurgarh - 124507, Haryana (INDIA). M.L.No.: MFG/MD/2021/000134
*Note: First digit in LOT (1 OR 2) refer to manufacturing unit.
E-mail : info@mediplusindia.com,
Web : www.mediplusindia.com

EU REP mdi Europa GmbH,
Langenhagener Str. 71, 30855
Hannover, Langenhagen, Germany
ingrid.sander@mdi-europa.com